Coding the meter

Code meter before using it for the first time and every time you change to another vial of OneTouch® FastTake® Test Strips.

STEP 1

Enter code mode.

Start with meter turned off. Insert OneTouch® FastTake® Test Strip to turn on meter. All display segments will appear followed by date and time. Code number will appear next.

STEP 2

Match code numbers.

Compare code number on meter display with code number on test strip vial. If they match, begin testing. If they do not match, follow Step 3.

STEP 3

Code meter.

Press C button until code numbers match. After  appears, OneTouch® FastTake® System is ready for testing.


Display symbols and messages

Blood glucose test result with suggestion to check ketone levels.

Average of last 14 days (64 tests performed).

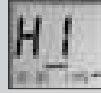Blood glucose level higher than 33.3 mmol/L (600 mg/dL).

Blood glucose level lower than 1.1 mmol/L (20 mg/dL).

Temperature too low for testing.

Temperature too high for testing.

Battery power too low for testing. Replace batteries immediately.

Note: Millimoles per litre (mmol/L) is the standard unit of measure in Canada. Select mmol/L whenever you reset the meter.

Problem with meter.*

Meter electronics problem.*

Blood sample applied before  appeared on display.*

Test strip damaged, removed during test, or inserted improperly.*

* Retest with a new test strip. If problem persists, call the OneTouch® Customer Care Line at 1 800 663-5521.

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Accurate results in just 15 seconds.

STEP 3



Apply sample. Touch blood drop to TOP WHITE EDGE of test strip. Do not push your finger against test strip. Hold blood drop to TOP WHITE EDGE of test strip until blood has completely filled confirmation window.

STEP 2

Insert test strip. Meter turns on automatically.

STEP 1

Note: With clean, dry hands, removing it from the vial or anywhere on its surface when you may touch the test strip inserting it into the meter.

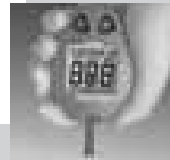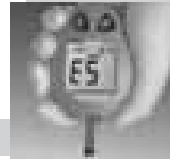Test summary

Exit memory mode. Press M button to turn off meter.

STEP 3

Recall test results. After three seconds, most recent test result with date and time appears. Press C button once and next most recent test result will appear. Continue pressing C button to recall last 150 test results in order.

STEP 2

Enter memory mode. Start with meter turned off. Press M button. When 14-day average appears, you are in memory mode.

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